

# COMMUNICATION CENTER 1-800-262-5221

AOA Wisconsin 13931 Spring Street Sturtevant, WI 53177 Fax 262-886-6879 AOA Connecticut 6 Niblick Road Enfield, CT 06082 Fax 860-741-7655 AOA California 341 E. First Street Calexico, CA 92231 Fax 760-357-9488

Dr. Pitts

INTERNATIONAL # 262-886-1050

### **BRACKET & TRAY SYSTEM STANDARDS**

Brackets are supplied by office, and sent to AOA Mesh Pad will have <u>Custom Composite Base</u>

### **TWO PIECE SYSTEM**

Clear soft Vinyl w/Clear Hard Support Tray

PL	AC	EM	EN	IT C	PT	10	NS
----	----	----	----	------	----	----	----

Standard	l	Dr.	Damor
----------	---	-----	-------

# PLEASE INDICATE WHERE TO SECTION

Do Not Sect	ion Tray
-------------	----------

- Section At The Midline
- Section Distal to Cuspids

## PLEASE WRITE SPECIAL INSTRUCTIONS

PLEASE SHIP EXTRA:
--------------------

- SHIPPING BOXES
- PRE-PAID BAGS
- PRESCRIPTION SHEETS


	ICE C	MLY
LAB	USP	


# Labial Indirect Bonding Rx

Dr	Acct #
Address	
City, State, Zip ————	
Patient	
	Fax #
E-Mail	
Shipped ———P	lacement Date
	E 1-2 DAYS BEFORE ACTUAL INSERTION DATE)
NOTES:	

Standard Slot Relationship Are Illustrated Below.
Please Indicate Any Changes On Chart
Circle Teech To Be Banded By Doctor. Place "X" On Teeth Not To Be Tx.



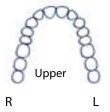

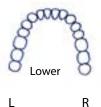

White and Yellow: Laboratory Copy Pink: Doctor's Copy

001-214-Fil REV D

Indiract	Pandina	with I	ight Cure	Adhacive
indirect	Bonaina	with i	iant Cure	Adnesive

Indi 1.	rect Bonding with Light Cure Adhesive Paint OrthoSolo® (Ormco) on the back of the brackets in the Indirect Bonding trays. Blow excess off after 30 seconds. Immediately place a small amount of Enlight LV® (Ormco) light cure adhesive on back of the brackets. Place under an orange cover or remove from light source.		
2.	Clean all teeth using pumice and prophy cup and rinse.	Labo	ratory Use Only
3.	Place dry angles and Nola lip retractors.		
4.	Micro-etch labial surfaces of posterior teeth.		
5. 6.	Etch labial surfaces of posterior teeth.	IB	
o. 7.	Rinse thoroughly. Dry teeth with tooth dryer.	ID	
, . 8.	Paint OrthoSolo® on teeth. Do Not Cure.	IB	
9.	Place Indirect bonding tray onto mandibular arch seating firmly.	ID	
10.	Cure each bracket at least 10 seconds at the Burst Cure setting with the Optilux 501® (Demetron Corp) while pressing moderately on each bracket.	IB	
	Repeat for maxillary arch.		
12.	To remove trays, first remove hard outer tray, then gently peel soft tray from each bracket while holding a finger on bracket. Repeat on maxillary arch.	IB	
13.	Remove excess adhesive with slow speed and check interproximal contacts with floss.		
		IB	
		IB	
		ID	
		IB	
		טו	
		IB	
		IB	
		IB	
		IB	
		IB	
		טו	
		SHIP	

RECEIVIN	NG	SHIPPING	
OPEN .		PULL	-
ID .		LAYOUT	
DATE .		PACK	
ENTER .		SHIP	
PULL .		CHECK	